



# MPath Spotlights – Image Browser

## Creating slide shows

Using Image Browser, you can document your work as you go, creating on-the-fly presentations to share or post to a website.

## Capture images

All MPath viewers include a Screen Capture option, which can be used to grab an image of whatever is currently displayed in the viewer.

To capture images:

1. In the viewer, from the File menu, choose Screen Capture, or click the Screen Capture tool button.
2. Enter a caption for the image and click OK.

Captured images are stored in a ScreenShots directory in the project directory.

## Organize the images

To organize the screen shots into a slide show, complete with captions:

1. In the main window, under Viewers, click Image Browser.
  - Image Browser lists all the images in the project directory, not just captures in the ScreenShots directory.
  - As you select images in the browser list, the image and its caption are displayed.
2. To change a caption, simply edit the caption. Use the View>Change Font option to change the appearance of the caption text.
3. Select images in the list and use the Edit menu options to Delete or Add images, and change the order of images (Move Up, Move Down) in the image list.
4. Use the View>Slideshow Mode option to only display the image; use the arrow buttons to change slides.

To save the current image list as is, use the File>Save option.

## Other things to try

You can export the slide show to a browser-viewable format:

1. From the File menu, choose Export Slideshow.
2. Enter a Directory name. The directory will be created in the SlideShows directory in the project directory when the slide show is exported.
3. Set the slide show options and click OK.

To display the slide show in any web browser, simply point your browser to the index.html file in the exported slide show directory.


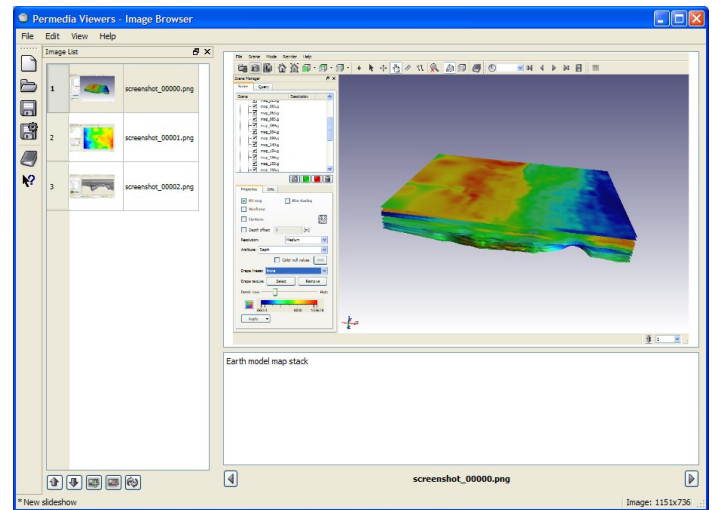


Image Browser

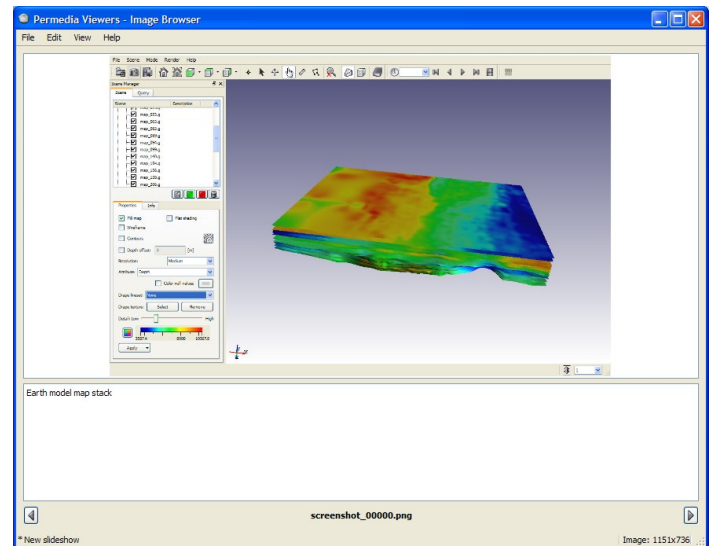

Slideshow mode